

MagniLink Student



MagniLink Student

More and more people with low vision are on the move and want to take their computers and video magnifiers with them. MagniLink Student is a portable video magnifier with combined reading, close up and distance viewing. It weighs only 1,5 kg (3,3 lbs) which makes MagniLink Student ideal for use in school and work environments where portability is important.

Ready, set, go

The unit is mounted in a couple of seconds and the connection to a computer is easily done through USB. MagniLink Student is powered from the computer and no external power supply or battery is necessary.

Ease of use

All functions are handled directly from the computer keyboard, or from a separate control box.


It is very easy to change between reading, close up and distance viewing. The camera head also rotates both vertically and horizontally, providing more flexibility for the user. The camera picture can be set to any preferred size and position and there are also preset split modes.

Student Professional & Student Addition

MagniLink Student is available both with a High Definition camera and a Standard camera.

The Student Professional comes with a High Definition camera which provides superior picture quality with crystal clear sharpness. The Student Professional provides a great overview with a magnification range of 1,5 - 55 x using a 15" monitor. With the great overview and the unique picture quality the Student Professional is ideal for use in school and work environments.

MagniLink Student Addition comes with the Standard camera and the possibility to connect directly to TV.

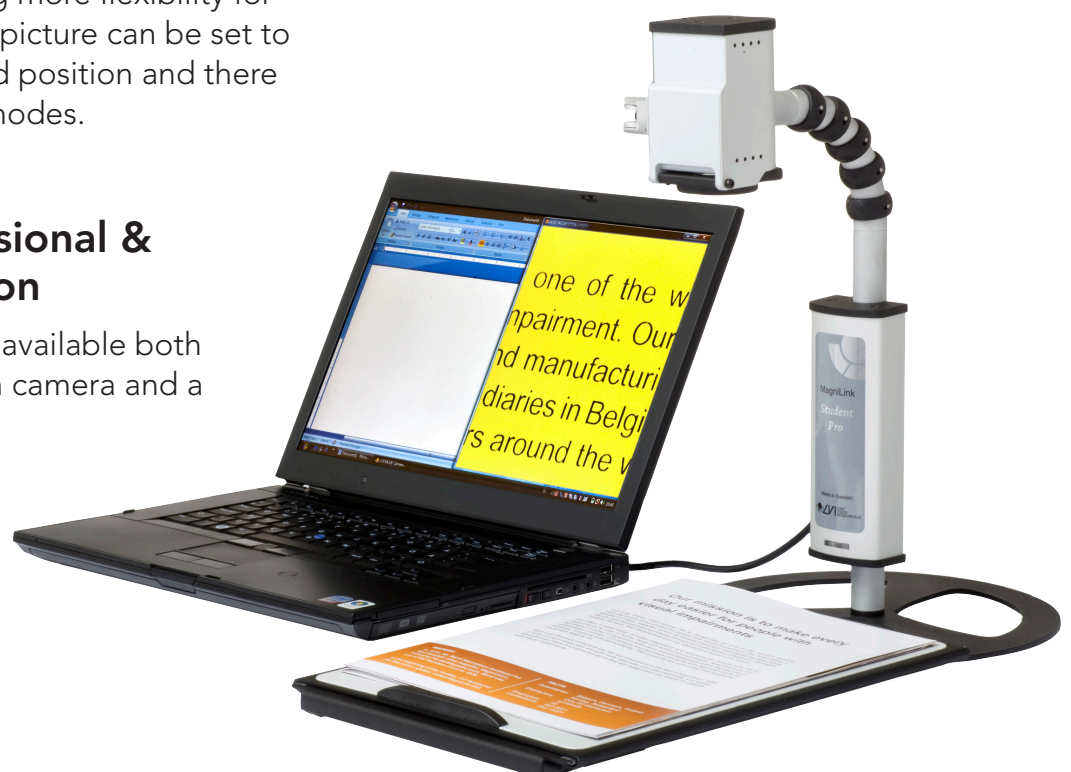

Features & benefits

- Superior picture quality thanks to the highest refresh rate on the market when used on a computer, allowing smooth reading.
- Modern design with a foldable and flexible arm, which makes the unit easy to carry and easy to place.
- The fast setup time allows the user to quickly mount the device when changing class- or meeting rooms.
- Memory settings for the different viewing modes, the unit will automatically adjust the mode and magnification to previous setting.
- The MagniLink S software is easily controlled by logical short commands, for instance simply press F5 to take a snapshot.

Further possibilities

- MagniLink Reading Table is a portable reading table for easy reading.
- MagniLink Docking Station is a high quality stationary reading table. Use the docking station at work and only bring the camera for meetings or visits to clients.
- MagniLink S Software can be upgraded when a new version is released. Visit www.lvi.se to learn about the latest features and upgrade your MagniLink Student.



TECHNICAL SPECIFICATIONS

Model	MLSTPRO or MLSTA
Connection	Laptop, desktop computer or external monitor MLSTA can also be connected to a TV
Camera	Standard: 640 x 480 High Definition: 1280 x 720 (720p)
Magnification	1,5–55 x (15,4" monitor)
Colour modes	Natural colours, 6 different artificial colours (positive/negative)
Weight	1,5 kg (3,3 lbs)
Refresh rate	PC: 60 frames/second External monitor: 60 frames/second
Compatibility	ZoomText, Lunar, LunarPlus and Supernova
Requirements*	Microsoft® Windows® XP with Service Pack 2 or Microsoft® Windows® Vista

* MagniLink Student has been tested on several computers without problems, but there is no guarantee that it will work on all computers following the specification. For recommendations contact, your local distributor or LVI Low Vision International AB.

Local distributor:

LVI Low Vision International
Verkstadsgatan 5
352 46 Växjö, SWEDEN

Tel: +46 (0)470 72 77 00
Fax: +46 (0)470 72 77 25
E-mail: info@lvi.se

WWW.LVI.SE